

9 Tube Professional Violet Ray System

Before Use

- Remove all metal jewellery (earrings , necklaces, bracelets etc)
- Ensure all items have been cleaned thoroughly before and after use.
- Items that touch the skin should be sterilized after each use.
- Glass tubes can be wiped clean with alcohol wipes.
- Always drink a large glass of purified WATER 10-20 minutes prior to any treatment. This is essential to hydrate the body.



Vacuum & Spray outputs

Violet Ray output

Galvanic outputs

Violet Ray Function

As mentioned it is essential to hydrate (drinking water 10 - 20 minutes before a session). Some people find (usually because they are dehydrated) that it can sometimes feel uncomfortable when they use the tubes on the skin!

Tubes

Connect the Violet Ray wand into the main unit carefully. Take care as the 5 pin plug only fits in one way. Do not force the plug in a way that would damage the 5 pins.



Clean and dry the area of the skin that is to be treated. Insert the desired tube carefully into the wand. Be careful not to bend or accidentally drop these glass tubes. They are very delicate and fragile. The tubes are vacuum filled with Argon or Neon gas that will become ineffective if broken from the metal bases.



1) Limb 2) Comb 3) Straight 4) Wart or point 5) Mushroom 6) Tongue 7) Roller 8) MWO or large area 9) Eye

General use of tubes – the tubes are not specific as they all function the same but the different shapes aid in using on certain areas of the body. The eye tube is useful for treating lines and wrinkles around the eye and has been reported to help vision (by Edgar Cayce). The comb tube is ideal for treating the scalp e.g. dandruff, dry scalp and restoration of hair follicles. The various other tubes can be used for lymphatic points, joints and/or arthritic points, massaging muscles, cleansing, toning and rejuvenating skin. You can also use on gums, nails and any points around the body. These are guidelines only, you can use the tubes on any area.

When using the tubes they do not have to touch the skin. There will be an electrical discharge or sparking effect when the skin is touched. This is in no way damaging. They can be used just above the skin if that is more comfortable.


Move the tube over the area you want to treat not staying in any one area for an extended period of time. As always, start with a shorter time period and adjust the intensity level to your own comfort levels.



After using clean the glass tube with an alcohol wipe and store the tube safely.

Vacuum & Spray Function

For the vacuum function – connect the red rubber tube to either one of the 3 glass tubes

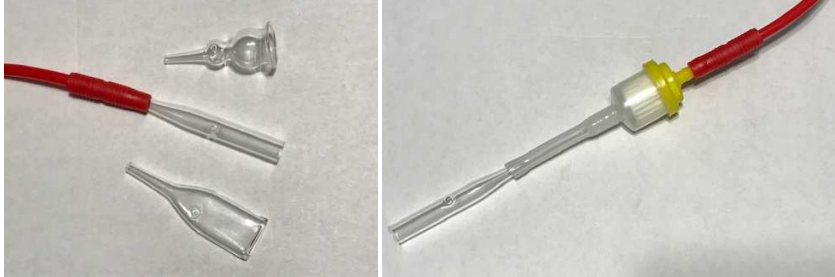


For treatment insert the glass tube directly into red tube. If desired you can insert the glass tube into the paper filter with yellow cap. Please note the paper filters are not available for replacement.

Insert the other end of the red tube into the main unit (left side of main unit).

Prepare the skin by using a natural oil (e.g. almond or olive oil) to allow smooth movement of the vacuum attachment. This helps with lymphatic drainage of the skin.

This picture on the right shows the starting points and direction of the vacuum function as you sweep the attachment across the face in the direction of the arrows

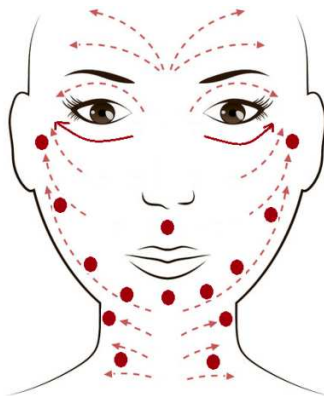

Adjust the vacuum intensity with the dial on the main unit and use your forefinger on the glass attachments to cover the air hole to begin the vacuum process.

For the Spray function – Connect the black tube to the spray bottle and the other end into the spray function on the main unit.



Fill the spray bottle with either distilled water or any scented water of your choice (e.g. rose water)
Switch on the Spray button on the main unit and hold your finger over the hole on the top of the spray bottle (as shown in the picture above on the right) to allow the spray to emit from the small hole on the side of the metal lid.

Galvanic Skin Option

Used for the toning of skin, the Galvanic option uses soft electric currents to stimulate, refresh and rejuvenate the skin.

Insert the electrode attachment into the far right output of the main unit. The recipient must hold this electrode in their hand while the roller or acupuncture attachment is being used. This completes an electrical circuit so the galvanic function can operate. Raise the intensity of the dial to suit the recipients comfort levels.



Acupuncture wand

Hand held electrode

Roller wand

ALWAYS clean all apparatus before and after each use. Clean the glass tubes and metal wands with alcohol wipes – DO NOT EMERSE THEM IN WATER!

Wash the rubber tubes by soaking them in warm/hot water with a sterilizing solution (e.g. baby bottle sterilizer) and hang to dry.

Pictures and video on the Violet Ray use can be found here:

<https://alteredstatesinstructions.com/9tubevioletray/index.htm>

Additional pictures and video will be added periodically