

# BioPhoton LS 5.0™

Added 3.5mm Input Port  
(not shown in all pictures)

## Full-Spectrum Photobiotic Nutritional Technology



**ADVANCED NOBLE GAS-TUBES  
With RADIO-WAVE SOUND PROBE  
HARMONIC FULL-SPECTRUM  
LIGHT, SOUND, OZONE &  
LYMPHATIC STIMULATION  
WITH ADDED INPUT PORT FOR  
ADDING CUSTOM PROGRAMS**

**COPYCLAIM © 1998-2021**

1998-2012 RESERVING ALL RIGHTS. 26 CFR (1939) 9.22(b)-1, UCC 1-308, UCC 1-103.6, 28 USC 1746(1), UCC 503  
NO HEALTH CLAIMS ARE MADE NOR IMPLIED

**The Lymphatic System** is a complex network of lymphoid organs, lymph nodes, lymph ducts, lymph tissues, lymph capillaries and lymph vessels that produce and transport lymph fluid from tissues to the circulatory system. The lymphatic system is a major component of the immune system. It has three interrelated functions: (1) removal of excess fluids from body tissues, (2) absorption of fatty acids and subsequent transport of fat, as chyle, to the circulatory system and, (3) production of immune cells such as lymphocytes (e.g. antibody producing plasma cells) and monocytes. Cells have an electromagnetic charge that either helps to bind them together or keep them apart. When they die or become damaged, their protein structures, many of which are damaged, must be removed by the lymphatic system. These nonfunctional proteins have electrical properties of attraction; they then bond to other proteins and water forming clusters that clog up the lymphatic system, potentially resulting in unbalanced conditions of abnormal growths or chronic inflammation - which in turn may lead to serious disease conditions in the body. **Cold gas light photons and extremely low-energy electromagnetic frequency patterns** help the cells to separate from each other and their accumulated fluids. This enables proper lymphatic drainage, allowing the body to rapidly rid itself of inflammation, swelling, abnormal growths, viruses, bacteria and fungus that feed on the protein wastes, food additives, and other chemicals, that attach to the protein waste.

Do You Suffer From TOXIC OVERLOAD? Symptoms may include: Headaches, Fatigue, Intestinal Problems, Back Discomfort, Depression, PMS, irritability, Inflammation/edema, etc.? All of these issues and many more may be often caused by toxic waste (i.e. chemicals, heavy metals, environmental pollutants, etc.) that builds up in the body due to of a congested or overburdened lymphatic system. A clogged or sluggish lymphatic system prevents the body from circulating vital fluids and eliminating toxic waste and dulls the immune system's response. This makes us vulnerable to swelling, infection, pain and a whole host of dis-eases. Unlike the blood system, the lymphatic system is not closed and has no central pump. Lymph movement occurs slowly with low pressure due to peristalsis, valves, and the milking action of skeletal muscles. Like veins, lymph travels through vessels in one way only, due to semilunar valves. This depends mainly on the movement of skeletal muscles to squeeze fluid through them, especially near the joints. Rhythmic contraction of the vessel walls through movements may also help draw fluid into the smallest lymphatic vessels, capillaries. Tight clothing can restrict this, thus reducing the removal of wastes and allowing them to accumulate. If tissue fluid builds up, the tissue will swell; this is called edema. As the circular path through the body's system continues, the fluid is then transported to progressively larger lymphatic vessels culminating in the right lymphatic duct (for lymph from the right upper body) and the thoracic duct (for the rest of the body); both ducts *drain* into the circulatory system at the right and left subclavian veins. The system collaborates with white blood cells in lymph nodes to protect the body from being infected by cancer cells, fungi, viruses or bacteria. This is known as a secondary circulatory system. Normally, the lymphatic system is activated by the movement of our body's magnetic field, deep breathing, exercise and activity, through that of the earth's magnetic field. It is essential to keep energy and fluids moving in the body so that its own natural intelligence may operate in full healing capacity. In addition, each cell must be enlivened with its own unique frequency, ideal energy-state, and be fully connected to the life-force energy throughout the rest of the body. Lymphatic system movement is essential in detoxifying the body, supporting the immune system and maintaining homeostasis.

Studies regarding the importance of the lymphatic system explain, most health practitioners' focus on feeding the body or enhancing its ability to eliminate, but the idea of transport, the process in between digestion and elimination, is many times forgotten. If nutrients can't get past the digestive tract and be transported to the area of need, then your nutritional supplements will not be absorbed properly.

**Pulsed Electrical Magnetic Frequency (PEMF) and Photon Technology** has been successfully utilized for over 50+ years to treat a wide variety of the major "dis-eases" that still plague the world today. Brilliant men like Nikoli Tesla, Royal Raymond Rife, Dr. Albert Abrahms, M.D., Georges Lackhovsky and Ed Skilling proved this type of technology very effective in the early 1900's. Stimulating the lymphatic system by the use of this technology has been widely recognized, used, and well documented, in Canada, Mexico, Europe, South America and Asia, and elsewhere, as an effective aid. PEMF and photon technology communicate with the body on a cellular and sub-cellular level to recharge and rejuvenate dysfunctional cells and body systems back to their optimal state. An electronic circuit that drives the energy by means of a pulse repetition rate in the audio spectrum produces the high voltage energy. This pulse rate has proven over a period of time to be of benefit to the human body as a form of healing energy much like normal exercise or deep breathing.

**Photons** are the light equivalent of electrons; ROS (Reactive Oxygen Species) electron scavengers. Without them electro-chemical reactions would not take place and life would be impossible. Electrons receive plenty of medical attention in the form of antioxidants to combat reactive oxygen species (free radicals). Photons may help restore the natural dispersion (via zeta potential) between molecules and facilitate decongestion of the lymphatic system, eliminate toxic waste from the body easier and more effectively and help to restore its ability to transport nutrients. Physics dictate **photons** are the driving source of our metabolism, whereas chemicals are only secondary stimulants and patterns of light information in this dimension. Therefore, ordinary nutrition through substance falls short in meeting our physical demands for full energy/information directly from the sun.

**Cold gas light photons** transfer energy frequency patterns to the targeted lymphatic regions. This puts the correct charge back into the cells and breaks up the congestion caused by clusters of water and pooled protein wastes. This helps our body to remove unwanted toxic wastes such as Steroids that mimic hormones and attach to proteins, unnatural additives from our food supply, etc. With proper lymphatic drainage, they can be rapidly flushed out of your system. **Bio-photon light** applications achieve its results by rebalancing the electromagnetic charge that either binds cells together or keeps them apart. You may find it more EFFECTIVE to use the bio-photon light application BEFORE any form of cleansing. Subsequent kill offs and growths of mutated Candida, fungi, bacteria, and viruses may INCREASE congestion if the waste cannot be excreted as fast as it accumulates. Dumping toxins out of the cells into already clogged lymph channels will equal a slow and challenging cleansing process. Lymphatic massage may have minimal benefits if the lymph is impacted with waste. Lymph drainage may occur by means of light utilizing an extremely low-current, **cold gas light photons** to transfer energy frequency patterns to cells in the targeted lymphatic regions. Frequency modulation is accomplished with light. **The high-frequency EM radio-wave** type transmissions deliver micro-current energy deep within the body's tissues, removing energetic blockages and restoring balance. By transmitting enough frequencies, every cell, of every organ, nerve, and tissue, can find its own unique frequency, and vibrate in resonance, settling in to a new, balanced dynamic harmony. Proper biological function returns, and balanced, vibrant health may be restored, once the cell is re-energized.

**The use of PEMF technology has been reported to help:**

Alkalize body fluids

Promote relaxation and clarity

Promote the body's synthesis and/or utilization of Melatonin

Promote health by stimulating all cells in the body Open and cleanse the lymphatic and circulatory systems.

Improve balance and motor functions

Move energetic fluids through all channels for nutrition transport, energy and intelligence

Eliminate blockages and barrier issues,

Neutralize negative "thought-forms" and old stagnant energy patterns

Increase blood flow and enhance blood's ability to carry oxygen

Increase circulation and the speed of metabolism

Reduce and eliminate toxins, stress, pressure and/or swelling from injury and/or blockage(s)

Therapists use less physical stress & improve results when coordinating sessions with this technology.

Cosmetic uses include toning and accelerating removal of cellulite with body wrap, etc

Dissolve lumps, clots, protein and mineral deposits, uric acid crystals, etc.,

Devitalize and destroy the structural integrity of pathogenic microorganisms and eliminate low frequency hosts and other harmful invaders, i.e.: viruses, bacteria, yeast, molds, fungi, parasites, worms, sarcoma and mycotoxins

Re-polarize cells and molecules

Enhance neurological processes

Energetic balancing through complete, direct holistic photobiotic nutrition

Increase energy levels

Promote bone and connective tissue regeneration, even where degenerative conditions exist

Re-establish each individual cell to normal frequency and energy-state

Stimulate and strengthen the lymphatic, circulatory and immune systems through "Re-Polarity"

Recharge life-force energy and facilitate higher states of consciousness

**RF Technologies™** is a research and development company specializing in Energy Science technology. For over 15 years, **RFT** has collaborated with world-class researchers, scientists and healers of many disciplines from various countries around the world to bring you the most effective, natural, nutritional supplements and modalities. To learn more about "living water" go to [www.BeyondWater.org](http://www.BeyondWater.org) "Adult Stem Cell enhancement" or "Telomeres", go to [www.CelluStem.com](http://www.CelluStem.com)

**The BioPhoton LS 5.0™ and Sound Probe™** is a combination of the brilliant healing technologies of light, sound & frequencies in one compact unit, using extremely low-current, utilizing noble cold gas ionization, tube technology, (**cold gas light photons**), and a radio-wave sound probe transmitter.

The Photon glass electrode provides frequencies of light in the Infrared and Ultra-violet light spectrums, providing bio-available micro-amperage output current for a gentle yet penetrating delivery of full-frequency "photon energy", and Ozone (O<sub>3</sub>) production.

Athletes who would like to maintain consistent peak performance, body workers and health practitioners who wish to complement their existing therapy sessions, or anyone wanting to improve their overall well-being, at any age at any level of health, may find this technology very beneficial. The **BioPhoton LS 5.0™ and Sound Probe™** is a universal tool for complete photobiotic nutrition ideally used in conjunction with a balanced natural health program. It may be used for both home and professional use as a research technology. It is safe, effective and easy to use for treating the body, air and drinking water.



**The Main Unit** in a compact, durable case measures 3.0"x 6.3"x 9.8".  
(4.0 shown here)



**The RF Sound Probe "flat-pack"** measures 3/4"x4". It is part of the complete circuitry of the BioPhoton LS 5.0, so it must be plugged into the main unit always. It can be used with or without the glass tubes. It is the Radio Frequency component of the BioPhoton that broadcasts frequencies via radio waves at distances up to 65' or more. Its deeply penetrating radio wave component contains a full-spectrum sweep of harmonic frequencies and scalar subtle energies containing billions of coherent harmonics and subtle scalar energy, through radionics. The Sound Probe may be placed directly on bare skin, which allows it to act as an additional "receiving" electrode, or placed in a small room. The closer the Sound Probe, the more powerful the signal is to the body. The greatest saturation of energy/information occurs when the probe is on the body. It has been reported to improve a more restful sleep and very effective for deep pain, structural & dental issues, etc.



**The Ground Plate** is for additional grounding purposes when needed. The **copper** Ground Plate will not have much sensation when touching the skin unless the tubes are also touching the skin. If both tubes are being used at the same time, the output energy will be increased.



**The 1 1/2Amp 250v fuse** is located on the lower left, backside of the unit, in a quick release fuse holder. This unit draws very little energy from the power source and the fuse may never have to be replaced. High voltage output without tubes plugged in, faulty transformer or defect in circuit board may cause fuse to blow.



**The round glass tube electrodes** with easy grip foam handles, measure approximately 6-7" in length. They come in several shapes and gas mixtures for efficient, various applications. They are evacuated of air and filled under vacuum with Inert or Noble Gas(s) including combinations of Argon and Neon gas. Mercury is NOT added to our glass tubes, adding Mercury is for atheistic purposes only and adds the undesired toxic frequency and element of Mercury. High voltage at low amperage (milliamps) ionizes the gases to a "Plasma State". The high voltage will overcome the insulation of the glass tubes and impinge upon the skin of the person using the equipment producing a pulsed magnetic energy that helps to stimulate the lymphatic system with each pulse. Ozone (O<sub>3</sub>) is produced when the glass tubes make contact with the skin. The Ozone offers oxygenation at the cellular level, delivering oxygen for detoxification. If the tubes are close together on any area of the body the energy will be at the maximum between the two gas-tubes. Multi-Wave Oscillator effects with the glass electrodes add hundreds of thousands of RF frequencies. This modality can be used in a variety of experimental electrode combinations with "sending" and "receiving" glass tube electrodes and the RF Sound Probe and Ground Plate. This BioPhoton is mono polar. This means that only one wire is connected to each glass gas-tube electrode. When the body is absorbing the high voltage energy, the return source to complete the circuit, is via Body Capacitance. The body completes the circuit; it is like an antenna system radiating this energy through the air back to the source, which is the circuit within the enclosure that generates the healing energy. The disadvantage, if any, is if another person is near by the person that is using the equipment, should touch the bare skin of the user, both parties may experience a mild shock. Keep the glass tube electrodes clean with warm water and a gentle soap, or spray with "Silver Shield" colloidal silver (available thru RF Technologies) and wipe with a small terry cloth towel. Be very careful while cleaning the glass electrodes, they are delicate and are breakable.

## Glass Tube Electrodes



The **Pyramid/Spiral** electrode is best used as a "sending" electrode. It gives the most surface area contact for increased lymphatic stimulation and drainage and it also produces the highest levels of ozone.



The **Pyramid Massage** electrode is great for lymphatic brushing techniques on the entire body. This is also great for grounding purposes.




The **Straight** electrode may be used for the neck, groin, limbs, prostate (perineum, base chakra), lower abdomen or low back i.e. liver/kidneys and for grounding.



The **Y** electrode may be used to glide up and down the spine or fitting around the "Cavity" areas i.e., ear, nose or corner of eyes, between the Y shape. It may also be used at the lower abdomen or lower back for grounding purposes.



The **L** electrode easily hooks and tucks inside the pants at the belt line, on the bare skin, for grounding purposes. This shape offers many uses that are left to your experimental ideas.

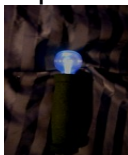

The **Regenerator** electrode is a small bulb shape that delivers a more powerful stimulating sensation, ideally for key lymphatic points. It is excellent for the neck, face and head as well as for epidermal toning. Great for running the Chakras, or Energy Centers of the body. Also known as the Endocrine System.



The All New "**Eye Tube**" electrode is similar to our Regenerator electrode with an extended neck for easier application. Also great for treating the male prostate externally.

## BioPhoton LS 5.0™ and Sound Probe™ Assembly Instructions



Connect Sound Probe into receptacle along right side, back of unit.



**Connect Glass Tube Electrodes:** Insert the male banana plug (@end of tube) into **red (positive output/sending)** and black (negative output/receiving) jacks on the front, lower left corner, marked "Tubes", on the main unit. Any electrode can be plugged into either the **red** or black jack for receiving and sending. The electrode plugged into the **black jack** is always used for "receiving" and the **red jack** is for **sending**. The two must be used together to complete the circuit. The jacks are "**High Voltage**" and are for the Glass Tube Electrodes **ONLY!** **Caution:** Do not plug the glass gas tube electrode banana plug into the jack if the **Neon Amber** lamp at center left is glowing or "On".



Connect the **copper sided Ground Plate black/male plug** into the black female banana jack on front; lower right corner, (under red button) marked "Ground". (For Ground Plate **ONLY!** **No Tubes should ever be plugged into this jack**)



Connect the **3.5mm jack to INPUT port** for adding custom programs and/or frequencies included with the 5.0 upgrade. 4.0 version can be upgraded to 5.0. ONLY use low level voltage inputs, such as a CD player, audio tuner, Light Box, etc..



Connect the power cord into transformer for power supply source.



Source: 120volt 60Hz external UL Listed transformer w/2-wire, 6' cord





**Connect power supply** into DC jack located behind unit, on right side.

Use surge protector in the "Off" position before plugging in the instrument to protect your unit from electrical surges. Plug the power supply cord into surge protector. You can use a 12V 2.1mm DC plug to the cigarette lighter port in most vehicles. Contact RF Tech to purchase this adapter.

## Placement Configurations

Place the **Sound Probe flat-pack** on the back just below the middle of the shoulder blades or place it on other sites (skin contact preferably) while using the glass electrodes in close proximity. You may also place the Sound Probe away from the electrodes such as the base of the spine, base of the skull, top of the head, between the knees (excellent), and all chakra points (especially Crown and Base). It is best to use the Sound Probe as the "ground-receiver", together with the sending and receiving glass tube electrodes at or near a trouble site.

**Note:** The Sound Probe may also be used alone (without using glass tube electrodes). When used alone for an overnight session, place the Sound Probe under your pillow or at the head of your bed.

Place the **Ground Plate, copper side touching skin**, at center below navel or under either bare foot. The **copper Ground Plate** is optional and does not have to be used with every treatment situation.

Place the **"receiving" glass tube electrode** in either hand, below the navel, lower back area, on the abdomen or the lower abdomen.

Place the **"sending" glass tube electrode** on Thymus, below the collarbones mid-sternum. **Always begin each session with 3-5 minutes of output power/sending**, moving the electrode from proximal to distal (from the center, out). This is to open up the main lymphatic thoracic duct into the main aorta so toxins can be quickly and easily flushed from the system. Be sure that electrodes always make contact with the skin.



The **sending** glass tube electrode is connected into the **red jack**. The **"receiving"** glass tube electrode is connected into the **black jack**.



Push the **red; main power, push on/off button**, on the upper right, front corner of the main unit. This powers the sound probe. The **Green Lamp** below the switch indicates that the power is "ON" and the electronic circuits are operating correctly and will blink at a 1 hertz pulse rate. If the **green** lamp ever stops its normal blinking or (pulsing), check to see if the power cord is plugged into an active wall receptacle. If the wall receptacle is active, then check to make sure that the fuse is good; an extra fuse is provided with every BioPhoton LS 5.0 and Sound Probe (BP-LS 4.0 version shown here)



Push the **black push on/off button**, on the upper left, front corner next. This powers the electrodes. The **Amber** (neon) lamp below this switch will glow in an on/off manner and the tubes will light with a color.



Move the front, center, Adjustable Intensity Dial clockwise and observe that both tubes are showing color and the intensity will feel stronger by increased levels of energy pulsation and ozone production. (BP-LS 4.0 version shown here)

Once all the surrounding lymph glands are opened, you may explore specific area(s) of concern by placing the "receiving" electrode near the site adjacent to the **"sending"** electrode or on the opposite side of the body. This way you pass the healing current through the site of interest. Regardless of where you localize the electrodes, the entire body is treated. When session is completed, set the tubes aside and return the Intensity knob to the "min." position and push the **"black push on/off button"** off first, and then push the **"red push on/off button"**, off.

**Group Usage:** Circle around, holding hands etc...Begin with someone holding the "Sending" electrode at one end, and someone else holding the "receiving" at the other end. Begin session with the suggested 3-5 minutes of **output power/sending**.

**Massage Therapist:** Use of the BioPhoton LS 5.0 is similar to a massage; it stimulates the increase of blood flow and relaxes muscles. Coordinating the BP-LS 5.0 prior to/with your massage therapy session is an ideal way to augment and intensify positive results. After 3-5 minutes **output power/sending**, start at soles of the feet and slowly stroke upward, toward crown chakra, then follow with massage therapy. Very rewarding!



**Always drink plenty of good living/structured water.** Remember the suggested water consumption is  $\frac{1}{2}$  ounce water, per 1-pound body weight. For example,  $\frac{1}{2}$  of 200 lbs.=100 oz of water each day. Your drinking water can be detoxified and ozonated by immersing the glass portion of straight electrodes into a container of water for 5 to 30 minutes. You can also remove magnetic imprinting on your water within 80-120 seconds, giving you what is known as "Blank Water" that is ready for reprogramming and information. "***Aquanutria Coral Calcium***" from Japan (*available thru RF Tech at [www.WaterIsLife.us](http://www.WaterIsLife.us)*) is a very effective aid in the detoxing process, cleaning and reprogramming/structuring of your drinking water. Our latest innovation is using BioCeramics to enhance your H<sub>2</sub>O experience. To learn more on what can be done using BioCeramics to heal your drinking and bathing water, go to: [www.BeyondWater.org](http://www.BeyondWater.org).

## CHOOSE BETTER HEALTH

THE PURPOSE OF THIS TECHNOLOGY IS TO IMPROVE HEALTH AND CREATE A SENSE OF WELLNESS. TO DO THIS, CHANGE HAS TO OCCUR, INCLUDING DETOXIFICATION AND REORDERING. THIS TECHNOLOGY WILL NOT ELIMINATE THE BODIES NEED FOR CLEAN ENVIRONMENT OR DETOXIFICATION. PLEASE KEEP THIS IN MIND AS RESEARCH PROGRESSES. THE INTELLIGENCE OF EACH PERSON'S BODY WILL USE THIS NOURISHING AND COHERENT ENERGY TO ITS HIGHEST GOOD AT ANY SPECIFIC POINT IN TIME. TO AN APPARENTLY WELL PERSON THIS WILL MEAN EVENTUAL DETOXIFICATION AND REORDERING INTO NEW LEVELS OF WELLNESS. THIS MAY YIELD PERIODS OF DISCOMFORT GENERALLY OR AT SPECIFIC POINTS OF WEAKNESS WHILE THE STORED STRESSES AND TOXINS ARE BEING RELEASED. REMEMBER THIS STRESS/TOXIN RELEASE PROCESS SEEMS TO BE AN UNAVOIDABLE ASPECT OF CREATED WELLNESS AT THIS TIME. ENERGY LEVELS WILL ALSO INCREASE. TO A SICK PERSON THIS PROCESS WILL MEAN STRENGTHENING AND STABILIZATION UNTIL THE BODY IS BALANCED AND STRONG ENOUGH TO DETOX AND RELEASE THE LONG STORED STRESSES AND TOXINS WHICH CAUSED THE SICKNESS TO BEGIN WITHIN.

THE BEST RESULTS ARE OBTAINED WITH THE BODIES' NATURAL HEALING MECHANISMS INTACT AND UNINTERRUPTED WITH. ACUTE ISSUES RESPOND MORE QUICKLY THAN CHRONIC, AS WITH OTHER HEALING SYSTEMS.

I HAVE REALIZED THAT THE BODY IS A RUDDER MEANT TO TELL US IF WE ARE LIVING ON COURSE WITH OUR SOUL'S PLAN OR OFF COURSE. THIS TECHNOLOGY WILL NOT CHANGE YOUR SOUL'S PLAN. IF YOU CONTINUE TO LIVE OFF COURSE THIS TECHNOLOGY WILL HELP YOU ONLY TEMPORALLY. HOWEVER, IF YOU WISELY CHOOSE TO TAKE FULL ADVANTAGE OF THIS TECHNOLOGY TO HELP YOU GET ON COURSE OR TO STAY ON COURSE, OF COURSE IT WILL TREMENDOUSLY HELP YOU TO REALIZE YOUR DIVINE DESIGN. SO USE THIS LOVING, NOURISHING, ENERGY WISELY AND ALLOW IT TO HELP YOU MANIFEST YOUR FULL POTENTIAL.

PLEASE BE FULLY AWARE OF THESE THINGS. THIS MACHINE IS NOT A "MAGIC BULLET" WHICH CAN SAVE YOU FROM YOURSELF. USED IN CONJUNCTION WITH CONSCIOUS AWARENESS, DIET AND SLOW DETOX, THIS NOURISHING ENERGY CAN BE AN INVALUABLE ASSET IN CREATING WELLNESS THROUGH CHANGE.

FOR WELLNESS IT IS IDEAL TO KEEP THREE THINGS IN YOUR AWARENESS:

1. THE PATTERN IN CONSCIOUSNESS THAT KEEP RE-CREATING THE DIS-EASE: THE ONLY WAY I KNOW OF TO CHANGE THIS FUNDAMENTAL CREATIVE PROCESS IS TO DO A TRANSCENDING MEDITATION, WHICH WILL DISSOLVE THE STRESS AT ITS CAUSE. IF THIS IS NOT DONE THE PATTERN WILL CONTINUE TO MANIFEST EITHER IN ITS ORIGINAL FORM OR A NEW GUISE.

2. **IN:** BE CONSCIOUS THAT YOU TAKE IN EVERYTHING THAT IS NOURISHING, LIFE SUPPORTING AND ALIVE I.E. HEALTHY FOOD PROGRAM, GOOD LIVING/STRUCTURED WATER, SUNLIGHT, CLEAN AIR. HERBS, HOMEOPATHICS, COLOR, SOUND, ETC...

3. **OUT:** ELIMINATE NEGATIVE THINKING AND TOXINS THROUGH MEDITATION, COLONICS, ENEMAS, LIVER FLUSH, FASTING, EXERCISE, SKIN BRUSHING, BODYWORK, ETC...

THIS MACHINE IS OFFERED TO YOU AS RESEARCH TECHNOLOGY ONLY. CONSIDER THIS YOUR PERFECT HEALTH RESEARCH PROGRAM. THIS IS NOT A SIMPLE CAUSE AND EFFECT TECHNOLOGY I.E. TOASTER, REFRIGERATOR, WASHER, WHERE YOU PUT IN DIRTY CLOTHES AND GET OUT CLEAN CLOTHES. IF HEALTH WERE THIS SIMPLE IT WOULD BE WONDERFUL. THERE IS NO GUARANTEE ON RESULTS. THERE IS A WARRANTY ON THE MACHINE ITSELF. IF THE MACHINE STOPS WORKING THROUGH NORMAL USE, IT WILL BE REPAIRED FREE OF CHARGE UP TO ONE YEAR.

CONGRATULATIONS ON YOUR CHOICE TO RESEARCH THE POSSIBILITY OF CREATING YOUR PERFECT HEALTH REGIME. THIS IS SOMETHING WE CAN DO ONLY FOR OURSELVES.

# RF TECHNOLOGIES

[www.BioPhotonLight.com](http://www.BioPhotonLight.com)

[support@biophotonlight.com](mailto:support@biophotonlight.com)

(619) 609.7109

## **BioPhoton LS 5.0™ Order Information**

BioPhoton LS - version 5.0 w/Sound Probe and Ground Plate includes two straight tubes and added INPUT port and insulated, hard shell carrying case.

Sales Price \$3,200.00\*

Special pricing for Distributorships available – email: [support@biophotonlight.com](mailto:support@biophotonlight.com)

### **Warranty**

1 year manufactures workmanship parts and labor warranty included.

Warranty does not apply to the glass tubes with electrodes nor the power transformer.

Additional 1 Year Warranty 99.00 (Covers main unit and Sound Probe)

### **Additional Glass Tube Electrodes**

Straight 110.00\*

Pyramid Spiral 165.00\*

Pyramid Massage 145.00\*

Y-Shaped 145.00\*

L-Shaped 135.00\*

Mini Gas-Bulb "Regenerator" 135.00\*

Eye-Tube 145.00\*

Custom Tubes 135.00\*-\$195.00+

Insulated Carry Case 115.00\*

6 pc. Tube Kit w/carrying Case 795.00\* (\$870 value) \*\*Special 695.00

1-"Pyramid Spiral", 1-"Pyramid Massage",

1-"Y", 1- "L", 1-"Eye-Tube" and 1-"Regenerator" bulb

\* Shipping, handling and insurance, not included.

\*\* Special Price when purchased with the BP-LS 4.0, or BP-LS 5.0

Order by phone: (619) 609.7109 or (213) 633.7225 ext. 1

Online: [www.BioPhotonLight.com](http://www.BioPhotonLight.com)

[support@biophotonlight.com](mailto:support@biophotonlight.com)

We accept Postal Money Orders, E-Check, Wire Transfer, PayPal, American Express, Visa, MasterCard and Discover Card.

## **NOTICE AND DISCLAIMER**

We are legally obliged to inform you of the following: Holistic products such as pulsating magnetic wave therapy devices are not recognized in the human and animal medical sciences and therefore cannot make claims regarding the areas of treatments. The listed areas of treatment and application are not yet proven and still debated in medical circles.

RF Technologies deeply respects your religious and Constitutional Right to "Freedom of Choice" in your own personal health care. In order to share this technology, we must comply with FDA regulations and market the BioPhoton LS 5.0 and Sound Probe as research and experimental technology and must be understood as such. This technology is used for your own experimental and research use only, regarding the potential of bioactive frequencies. The private reports of users of these research instruments can be provided to you through your own research. Each person who purchases this experimental research technology from RFT is considered a "Researcher" and must be self-conscious/self-responsible in all possible ways.

RF Technologies assumes NO LIABILITY for the use and/or misuse of the BioPhoton LS 5.0 and Sound Probe. RFT, RFT employees, associates or distributors, cannot nor will not be held responsible/liable for the results achieved or not achieved from the use or misuse of our products or suggestions. It is not meant to treat, cure, prevent or diagnose any health condition(s), or to replace the advice of a health care professional, preferably one who is familiar with Natural Health & Alternative Therapy. No health claims are made nor implied. The information provided in this booklet and/or any of our associated websites, is for informational and educational purposes only. Statements made here have not been evaluated or approved by the U.S. Food and Drug Administration or any other federal regulatory agency.

RF Technologies reserves the RIGHT to present ANY information under the BENEFIT and AUTHORITY of the 1st Amendment to the Constitution of the united States of America and "SAID" AUTHORITY pertains to that Freedom of Speech and Freedom of the Press, ordained and preserved to the BENEFIT of ALL peoples. Also, this company reserves the RIGHT to present ANY Information under the BENEFIT and AUTHORITY of the \* 9th Amendment to the Constitution of the united States of America to provide TRUTHFUL and ACCURATE Information about our "Holistic Devices" WITHOUT receiving PRIOR approval from ANY governmental agency or agencies (i.e. - F.D.A., F.T.C.). This RIGHT is being FREELY used NOW and will be used AGAIN

This consent and waiver agreement is effective on the date shown below by and between RF Technologies (the "SELLER") and the undersigned buyer of goods and/or services (the "BUYER"). BUYER acknowledges, understands, and agrees that:

1. The products, equipment, and devices (the "Products") sold by SELLER and purchased by BUYER are for cosmetic use, are not approved by the U.S. Food and Drug Administration for medical treatments, and are not intended to be used for medical treatment.

2. SELLER is not a licensed physician. Any use of the Products, consultation, or other services provided by the SELLER to the BUYER are for demonstration purposes only, and are not intended to be nor may they be construed as a medical treatment or a substitution for medical treatment by a licensed physician.

3. BUYER has consulted with a licensed physician regarding any and all medical conditions that he or she may have, or has independently and voluntarily, without compulsion, coercion, or compensation failed to have such a consultation.

4. BUYER voluntarily agrees and consents to participation in any and all demonstrations of the Products requested by the BUYER and administered by the SELLER. Any demonstration shall be terminated immediately upon the request of the BUYER.

5. SELLER'S WARRANTY IS STRICTLY LIMITED TO THE WARRANTY PROVIDED BY THE MANUFACTURER OF THE PRODUCT, AND NO OTHER WARRANTIES ARE PROVIDED, EITHER EXPRESSED OR IMPLIED, FOR PURCHASE, USE OR DEMONSTRATION OF THE PRODUCTS. SELLER SHALL NOT BE RESPONSIBLE FOR, AND BUYER WAIVES, ANY CLAIM, CAUSE OF ACTION, SUIT, OR OTHER REMEDY AGAINST SELLER OR LIABILITY THEREFORE EXCEPT AS SUCH CLAIMS, CAUSES OF ACTION, SUITS, AND OTHER REMEDIES OCCUR AS A RESULT OF SELLER'S SOLE NEGLIGENCE.

---

BUYER Signature \_\_\_\_\_

Print Name \_\_\_\_\_ Date: \_\_\_\_\_