


Multi-Pulsar



By

Altered States Ltd

www.altered-states.net

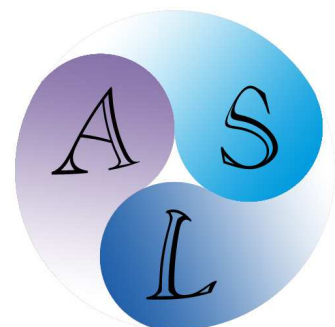

Congratulations and Thank You!

You have purchased one of the most powerful and advanced Magnetic Pulsar's available on the market today.

The **Multi-Pulsar** is designed to give you many years of reliable service.

Please read this instruction manual carefully before use

As with all electrical products **DO NOT USE:**


- *Near water
- *When driving a car
- *when operating machinery or heavy equipment
- *around unattended children
- *around metal objects

The Multi-Pulsar contains a powerful magnetic coil in the hand held paddle, so please keep away from homeopathic remedies, credit cards, security cards, audio and video tapes, computer discs, answering machines, computers and any other metal object that could be effected by a magnetic field.

The Multi-Pulsar is a very easy unit to use. It features bonus functions such as a Colloidal Silver Generator and a Bob Beck style Zapper along with the Magnetic Pulsar function. It can be operated on 110 volt and 240 volt power systems.

Best wishes on your research
from the team at Altered States

Made and Designed in New Zealand

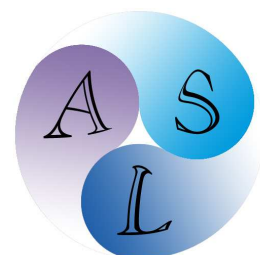

INDEX

Introduction.....	2
Index.....	3
Disclaimer.....	3
Technical Specifications.....	4
Description.....	5
Magnetic Pulsar.....	6
Zapper.....	9
Colloidal Silver.....	12
Troubleshooting.....	14

Disclaimer

Please note that our legal advisors have instructed us to issue the following disclaimer so as to comply with international trading agreements regarding this electronic frequency generating device.

Altered States products are sold for learning, self improvement, research and simple relaxation. No statement contained in this publication, and no information provided by any Altered States employee should be construed as a claim or representation that these products are intended for the use in diagnosis, cure, mitigation, treatment or prevention of disease or any other medical condition. The information contained in this publication is deemed to be based on reliable and authoritative reports. However, certain persons considered experts may disagree with one or more of the statements contained here. Altered States assumes no liability or risk involved with the use of the products described here. We make no warranty, expressed or implied, other than that the material conforms to applicable standard specifications. The publisher does not accept any responsibility for the accuracy of the information or the consequences arising from the application, use or misuse of any of the information contained herein, including any injury and/or damage to any person or property as a matter of product liability, negligence or otherwise. **No warranty, expressed or implied is made in regard to the contents of this material.** No claims or endorsements are made for any drugs or compounds currently marketed or in investigative use. **This material is not intended as a guide to self-medication.** The reader is advised to discuss the information provided with an authorized healthcare practitioner and to check product information (including package inserts) regarding dosage, precautions, warnings, interactions and contraindications before administering any drug, herb, electronic tool or supplement discussed herein.

Technical Specifications

The unit is based on a microprocessor controlled system which allows for multi function use, precise timing and efficiency.

It is compact with digital and multi coloured LED's which are easy to use and care for.

Power Supply: 36V 1A dc

Timer: Zapper, Colloidal Silver: 5-120 minutes (5 minute steps)

Pulsar, Fast Pulsar: 5-30 minutes (5 minute steps)

Intensity: Zapper: 5v-35v 4Hz (50% during AC square wave)

Colloidal Silver: 27v

Pulsar: Output 1 magnetic level only 1-75bpm 6500-7000gauss

Care of Device

- Wipe the Paddle with a dry non-abrasive cloth. DO NOT USE SOLVENTS OR LIQUID CLEANERS
- Clean Silver rods with a nylon scouring pad (similar to the one provided) or paper towels. It is normal for one of the rods to blacken during the process, this is silver oxide which should be removed.
- Clean the silver electrodes the same as above.
- Do not twist or tightly wrap the electrode cables, this can (and eventually will) damage the electrical cable.
- Keep the Silver rods and electrodes away from other metals, this is pure silver, so do not clean with metal scouring pads.
- DO NOT USE THIS DEVICE IF ANY DAMAGE IS VISABLE ON THE UNIT, Contact us immediately.

ALLOW MAGNETIC PADDLE TO COOL DOWN COMPLETELY BEFORE BEGINNING ANOTHER SESSION!

PLUG THE POWER CABLE INTO THE UNIT AND THEN INTO THE MAINS POWER BEFORE SWITCHING THE UNIT ON!

Description

The **Multi-Pulsar** is a consumer product designed to complement a healthy lifestyle.

Please understand that results will vary with individual users. Healthy lifestyle choices are important. Please read this manual carefully before using the **Multi-Pulsar**.

The **Multi-Pulsar** has 3 main functions:

- 1) Magnetic Pulsar
- 2) Colloidal Silver Generator
- 3) Zapper

Contents of Package:

- 1) Multi-Pulsar main unit
- 2) Power Adaptor (Multi voltage adaptor)
- 3) Product instructions
- 4) Lead with alligator clips (for colloidal silver function)
- 5) 99.99% pure silver rods (two)
- 6) Scouring Pad
- 7) Cable for small round gel pads (zapper electrodes)
- 8) 4 small round gel pads for zapper function

WARNING!

Before using any frequency generator, it is very important to drink water to hydrate your body. Ozonated water is particularly helpful.

DO NOT USE IF:

1. You are pregnant or trying to get pregnant
2. You use a pacemaker
3. You have metal implants

Magnetic Pulsar

ALWAYS plug power cord into unit and mains power **BEFORE** switching unit on.

DO NOT USE A POWER SURGE PROTECTOR ON THIS DEVICE!

Press the ON/OFF/SELECT button to switch on.

The display will show a hyphen (or minus) sign “-” on the left side of the display



Press the UP/DOWN button so the “-” is next to the PS (pulsar) and press ON/OFF/SELECT to program this.

The timer will flash. Press the UP/DOWN button to select the desired time 5,10,15,20,25 or 30 minutes.

Press the ON/OFF/SELECT button to program the timer selection.



Select the pulse rate (pulses per minute).

A “P” will appear indicating pulse selection.

Press the UP/DOWN buttons to select between 1-75 pulses per minute.

Press the ON/OFF/SELECT button to program the pulse rate.



The timer will begin to count down with a decimal point flashing indicating the timer counting down. The unit will switch off after the timer has finished.

The magnetic pulse will be emitted each time the counter on the left reaches “0”. It will count down from 9-0 and then pulse each time it reaches “0”.

The hand paddle contains the magnetic coil through which the magnetic pulse is released. With normal use the hand paddle may become warm. If it is uncomfortable turn the unit off and allow it to cool down before further use. DO NOT USE THE UNIT FOR MORE THAN 30 minutes FOR EACH SESSION? LET THE UNIT COOL BEFORE BEGINNING ANOTHER SESSION (approx 20-30 minutes).

The head of the paddle should be located on or as close to the body as possible. It does not need to touch the body. Clothing will not interfere with the magnetic pulse.

REMOVE ANY METAL OBJECTS FROM THE BODY BEFORE USING e.g. ear rings, belt buckles, jewellery etc or any metal that may interfere with magnetic frequencies.

Most people will not feel any sensation while using the Magnetic Pulsar function. This **does not** mean the unit is not working.

Magnetic Pulsar Continued

Session Length: Daily use varies from 5 minutes to 30 minutes. It is best to increase usage time slowly.

Length of Use: PMF's (Pulsed Magnetic Field's) are most effective if used daily or twice daily. Daily use is more important for results than using occasionally for a longer period. It is best to gradually increase the length of time the unit is applied.

Placement: The hand paddle may be kept on the one location one or more of the 30 minute sessions. It may also be moved at any time to a different location.

Water: It is very important to drink water prior, during and after using the unit. You must keep hydrated. Ozonated water is very helpful for this.

Stop Gradually: Reduce the amount of time each day for the last week of use. If using the suggested 4 part protocol it is suggested that you only stop one part at a time.

Maintenance Schedule: Here are suggestions for using the Magnetic Pulsar as a maintenance schedule:

- one or two sessions daily for two days per week
- one or two sessions daily for one week
- one or two sessions daily for one month repeated twice yearly.

Cautions: It is possible that using Pulsed Magnetic Field's can mask pain. It is wise to consult your health practitioners to determine the cause of the pain.

ALLOW MAGNETIC PADDLE TO COOL DOWN COMPLETELY BEFORE BEGINNING ANOTHER SESSION!

North/South

North polarity has an anti-clockwise rotation (negative polarity) giving a more 'normalized' effect on the body. We suggest the majority of Magnetic Pulsar sessions uses the North polarity. It promotes an alkaline environment.

South polarity has a clockwise rotation (positive polarity) giving more of a 'stimulating' effect on the body. It promotes an acidic environment.

Please research further benefits of North/South polarity.

Zapper Function

This is a very useful and well-researched program. Many people have reported beneficial results using this program as part of a detoxification protocol. Individual results may vary.

The Zapper runs 4 Hz equal pulse width bi-polarity square waveform to output. Intensity can be adjusted while the program is running. Timer should be selected from 5 to 60 minutes. The screen will display the timer counting down while a program is operating.

This information here is based on Bob Beck's original protocol for using the blood electrifier and colloidal silver generator. Please read this protocol to understand this procedure better.

According to his research, he found a way to kill most known viruses with a four pronged treatment: 1) Ozonated water 2) Colloidal silver water 3) blood purifier (zapping device) 4) Magnetic pulsar.

This unit is made to the exact specifications designed by Bob Beck. It can be used in for numbers #2, #3 as listed above. As all people are different, results may vary. These units are most effective when used according to BBs' protocol. Nothing is guaranteed. Put simply, the closer you match this protocol, the better results you may achieve.

Research from Harvard, MIT and Albert Einstein College of Medicine has shown that micro currents are known to eliminate viruses, parasites, fungi, bacteria and pathogens in the blood. This unit emits those same micro currents.

It has been a long known fact, that silver colloids greatly assist in eliminating known pathogens and prevents opportunistic infections.

Before zapping, it is very important to eliminate ALL herbs (including garlic), dietary supplements, alcohol, recreational drugs and medications. Stop taking these at least 24 hours prior to using the zapper (preferably 4-7 days before zapping). You may then take medicine or vitamins straight AFTER zapping. This is because zapping can allow medicines, drugs and herbs to enter the body's cells and it can actually magnify their effect 20-30 times so while you are zapping with these in your body, the supplements are amplified dramatically.

Zapper Electrodes

It is very important to drink at least 3-4 glasses of filtered or purified water (even distilled water) during the zapping process. This will help to flush your system and allow “detoxification” to occur. **YOU MUST DRINK PLENTY OF WATER WHILE ZAPPING** or you will not detoxify your body correctly.

It is important to place the electrodes correctly otherwise zapping will not be effective. Keep the space between the electrodes dry. Connect the electrodes to your wrist as specified, plug lead into ‘output’ **BEFORE** switching unit on. The intensity will always be set to **ZERO** when you begin (for safety reasons).

Everybody will have a different skin resistance, so people will ‘feel’ different effects. The different ‘intensity’ settings will accommodate most people. **BE AWARE** – you may not always ‘feel’ the zapping effect. This does not matter, the frequencies **WILL** be received, it’s just that your body will have adjusted its’ own skin resistance as a reaction to the zapping.

Bob Beck electrode kit

Contains 4 x reusable gel pad electrodes with connecting cable.

These are to be placed directly over the radial and ulnar arteries on the wrist.

After use, place the gel pads back onto the clear plastic pad they came on so they can be protected for further use.



Operating Zapper:


Connect main power cable to mains power and connect the zapper electrode to the OUTPUT on the front of the unit.

Switch power on by pressing the ON/OFF/SELECT

button

Press UP/DOWN to select ZP when the “-” is next to “ZP”.

Press ON/OFF/SELECT to program the zapper.

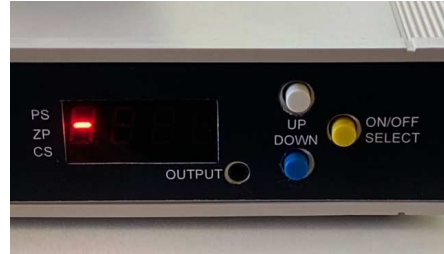

Press UP/DOWN buttons to select desired time 5,10,15,20,25 etc.

Press ON/OFF/SELECT button to program the time.



Next, select the Frequency. An “F” will appear on the screen. 4Hz, 10Hz, 20Hz,30Hz,40Hz press UP/DOWN/SELECT to program the desired frequency.



The zapper will begin with the timer counting down with a decimal point flashing indicating the timer operating.

The zapper will always start with the intensity at “0” (for safety reasons).

Press the UP button to increase the intensity to a comfortable level.

The intensity levels are 0-9.

The unit will switch off after the timer has counted down.

Colloidal Silver

Much has been written about Colloidal Silver and it's benefits. This is an important part of Bob Becks' suggested protocol.

The Colloidal Silver function has a fixed 27V output single polar (centre pin +) voltage and timer can be selected from 5 minutes up to 120 minutes. No features can be adjusted while this program is running.


Colloidal Silver water is very susceptible to light and temperature change, so it must be kept in dark amber bottles or opaque containers. A stable room temperature environment is ideal. Do NOT use plastic bottles (this will neutralize the colloids).

The silver rods supplied are 99.99% pure silver. A higher grade of silver of 99.99% purity can be used. Lower grade silver or sterling silver will contain nickel which is toxic and must NOT be used.

Connect the silver rods to the alligator clips; plug this lead to the “**output**” on the main device. Place the silver rods over the top lip of the glass (see picture below).

Use 250ml distilled water. ONLY USE DISTILLED WATER!

Place the silver rods 4-5cm apart. Keep them parallel and do not let them touch (otherwise it is likely to “short” the unit out and render it ineffective).

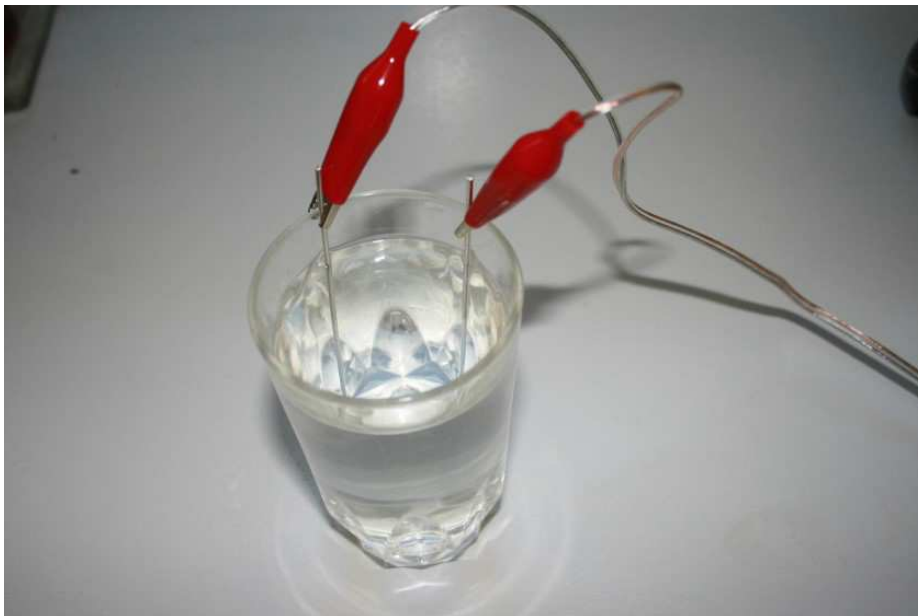

Instructions cont.


Connect main power cable to mains power

Switch power on by pressing the ON/OFF/SELECT button.

Press the UP/DOWN button to select

CS (colloidal silver)


Press ON/OFF/SELECT to program this function.

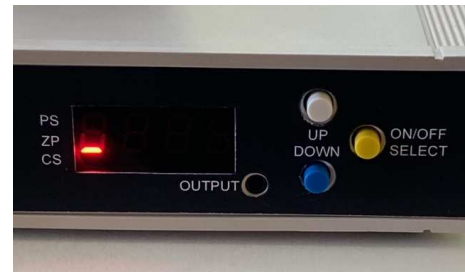

Press UP/DOWN buttons to select desired

generating time 5-120 minutes

press UP/DOWN/SELECT button to program the timer.

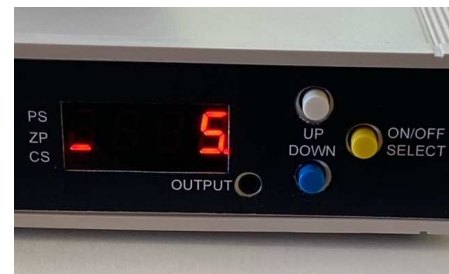

The unit will begin with a decimal point flashing to indicate the unit is operating.

The unit will switch off after the timer has counted down.

Remove the silver rods from the glass and clean with the nylon scouring pad.

Troubleshooting

If there is no power to the unit, please check the following:

- Ensure the power adaptor is fully plugged into the main wall power socket
- Ensure the power lead is firmly plugged into the rear of the Multi-Pulsar unit
- Try another wall socket

If you cannot feel the unit operating in Zapper mode:

- Ensure the Zapper electrode lead is plugged into the OUTPUT firmly on the front panel of the Multi-Pulsar
- Plug the Test Plug into the OUTPUT socket, turn up the intensity. The LED on the Test Plug will flash at the same frequency as the output frequency.

If the unit appears not to function correctly, please contact us immediately.

0064-9-828 9800 or altered@alteredstates.info

www.altered-states.net